


# CODAN Check Valves



The versatile programme  
of check valves



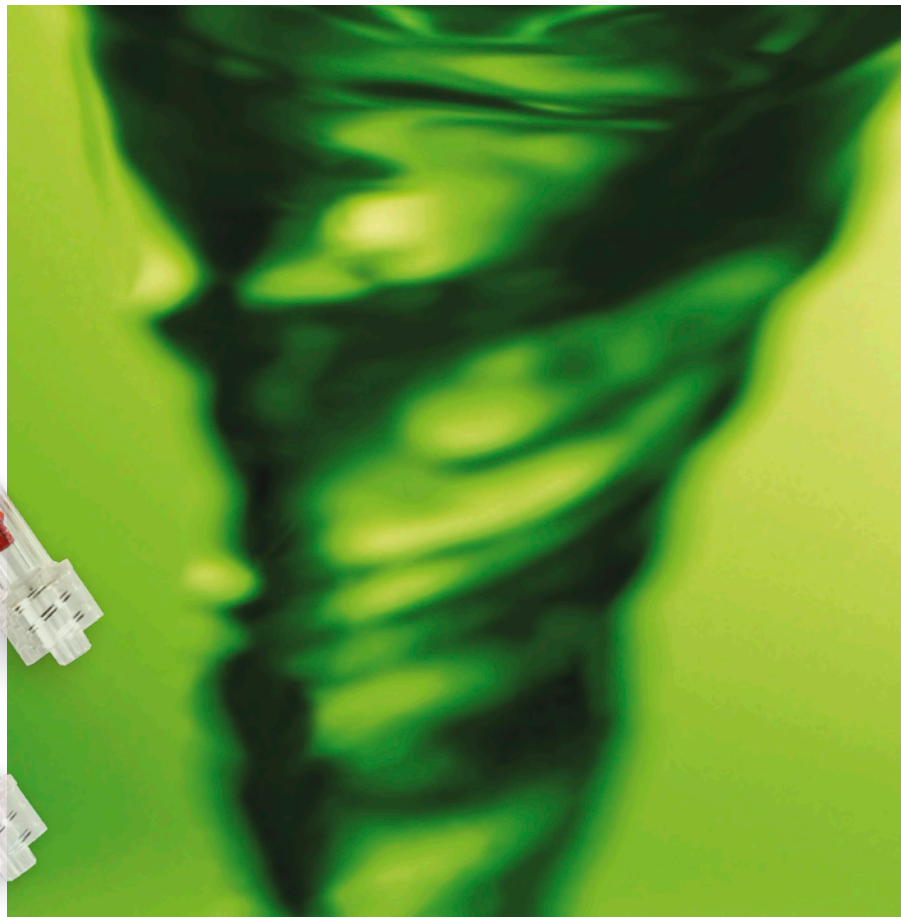

# CODAN Check Valves

## **Integrated safety**

CODAN Check Valves are a small detail with maximum importance for patient safety. As a firm or optional integral part in various CODAN products, the check valves reliably prevent reflux or back flow of infusion solution without additional measures and also ensure that blood does not flow back from the patient into the infusion system.

## **Small, perfect and safe**



### **Top-class engineering in a mini-format**

As with all CODAN products, the CODAN Check Valves also manifest decades of experience in research and development.

Precision in production ensures maximum reliability and ease of movement of the CODAN Check Valve. The valve opens upon the slightest conveying pressure. It closes automatically as soon as a counter-pressure results or the conveying pressure has dropped again. Even highly viscous infusions pass through the valves securely and reliably.

The CODAN Check Valves have a particularly low priming volume of 0.15 ml and can be integrated into all infusion systems without any problems.

Be it in standard infusion therapy on wards, in the operating theatre, in intensive care or in anaesthetics - CODAN Check Valves can be used in all areas of everyday hospital work.

# Versatility with precision

## For use in many versions

The matching solution can be found for each application. The CODAN Check Valve can be integrated into an existing system as an individual component. In addition, a variety of connection sets, extensions, connectors, completely pre-fabricated systems and special systems with CODAN Check Valves are available. All variants are available with differing tubing diameters. Connectors with CODAN Check Valves are available as double, triple or quadruple connectors.

## Connections

CODAN Check Valves in CODAN infusion set connections make handling particularly safe and simple - the less individual components are used, the more quickly and more reliably the infusion can be prepared and administered. Special importance is given to the Y-shaped injection site in safe and low-risk administration of cytostatics.

## Connection variants

The connections provided with a rotating collar which enable secure closure of the Male Luer-Lock without tube torsion are particularly comfortable and user-friendly. Of course CODAN also provides all standard infusion sets with integrated CODAN Check Valves.

R-LOCK,  
Female/Male  
Luer-Lock adapter  
with CODAN  
Check Valve

Inline adapter  
with CODAN  
Check Valve

Y-injection site  
with CODAN  
Check Valve

Double Y-site  
with CODAN  
Check Valves

Double connector  
with CODAN  
Check Valves

Rotating  
Male Luer-Lock  
adapter with CODAN  
Check Valve

Male Luer-Lock  
adapter with CODAN  
Check Valve

**CODAN**

# CODAN Worldwide



CODAN is known internationally as a manufacturer and supplier of medical transfer systems. The CODAN Companies have more than 1500 employees around the world.

The name CODAN is synonymous with reliability, quality and precision based on the know-how and experience gained from more than 60 years of research and development. Company-owned production facilities and sales companies around the world are a guarantee for efficient production, a tight-knit sales network and a first-class service for customers in the healthcare sector.

## CODAN Companies

CODAN Medizinische Geräte GmbH · Deutschland  
CODAN pvb Critical Care GmbH · Deutschland  
CODAN pvb Medical GmbH · Deutschland  
CODAN 11, S.A. · Portugal  
CODAN US Corporation · California · USA  
CODAN Inc. · California · USA  
CODAN NORGE AS · Norge  
CODAN TRIPLUS AB · Sverige  
CODAN Limited · Great Britain  
CODAN FRANCE Sarl · France  
CODAN Medical AG · Schweiz  
CODAN ARGUS AG · Schweiz  
CODAN BV · Nederland  
CODAN s.r.l. · Italia  
CODAN Medical GmbH · Österreich  
CODAN Medical ApS · Danmark  
CODAN DEHA ApS · Danmark  
CODAN MEDITECH s.r.o. · Česká republika

### CODAN Product range

- Infusion sets
- Transfusion sets
- Extension lines and manifold connectors
- Infusion and transfusion accessories
- Infusion filters and filter systems
- Neonatology/Paediatric products
- Withdrawal, preparation and administration systems
- CODAN CYTO®
- Chemoprotect® products
- Single use syringes
- Invasive blood pressure monitoring systems
- Infusion pumps
- Other CODAN Products

Compliance of the established quality management systems with the provisions of EN ISO 13485, the Council Directive 93/42/EEC and/or Regulation (EU) 2017/745 has been certified by the relevant, competent notified bodies:

### TÜV SÜD Product Service GmbH

CODAN Medizinische Geräte GmbH  
23738 Lensahn, Germany

CODAN pvb Critical Care GmbH  
85661 Forstinning, Germany

CODAN US Corporation · Santa Ana, CA 92704, USA

CODAN ARGUS AG · 6340 Baar, Switzerland

### Presafe Denmark A/S

CODAN Medical ApS · 4970 Rødby, Denmark

### Manufacturer

CODAN Medizinische Geräte GmbH  
Stig Husted-Andersen Straße 11 · 23738 Lensahn, Germany  
Tel. +49 (0) 43 63 · 51 11 · Fax +49 (0) 43 63 · 51 12 14  
codan@codan.de · www.codancompanies.com



**The decisive connection**